




ProMOTION[®]
EV1

QUICK START GUIDE

ProMOTION DEVICE

1 Battery Charging


Ensure your **ProMOTION EV1** battery is fully charged before use. Charge the battery using the supplied charger (fig 1) - insert the charger into the DC jack on the battery pack (fig 2)

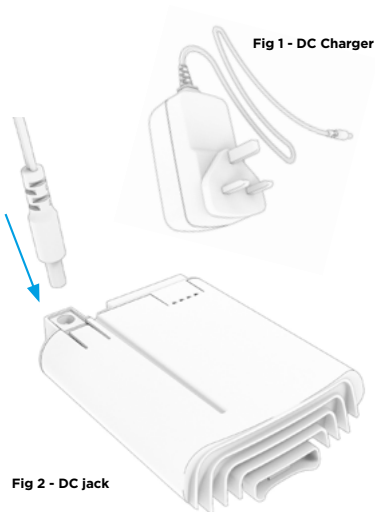


ProMOTION DEVICE

2 Assembly

Assemble the battery pack into the Therapy Unit.

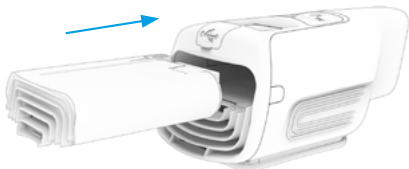


If using temperature therapy only, you can also assemble the Therapy Unit into it's securing strap for holding the device against the body

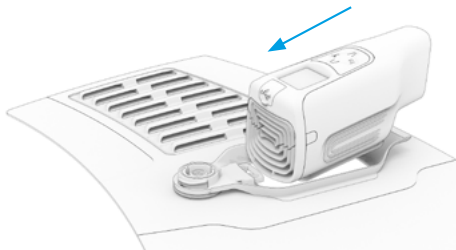


ProMOTION DEVICE

3 Assembly

If you're using compression therapy too, assemble the Therapy Unit and Pump Unit into the cradle.

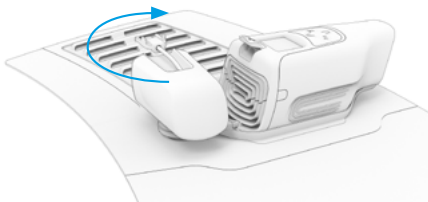


Connect the Pump Unit to the Therapy Unit with its USB cable.

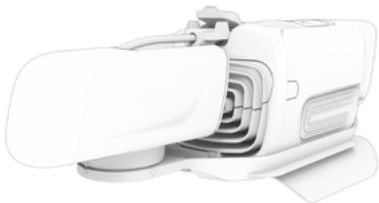

ProMOTION DEVICE

4 Location

The Compression Wrap should be located on the body. Align the Therapy Unit with the site of injury, and pass the wrap around the limb. The wrap should be secure enough not to slip, but not so tight that it is difficult to hook the handle in place.

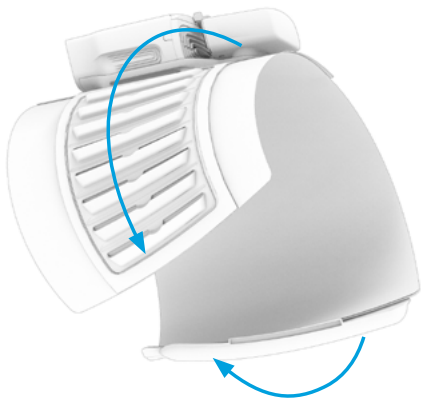

ProMOTION DEVICE


4 Navigation

Press and hold the power button until the device beeps, then release.

Selecting 'New' allows you to tailor individual settings for temperature, duration and compression levels.

Pressing Select allows you to progress to the next menu. Pressing 'Go back' returns you to the previous menu.


Press and hold the button to power the device off.

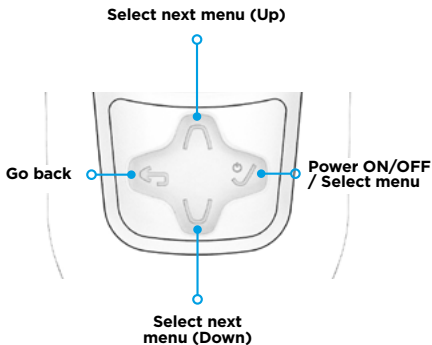

ProMOTION DEVICE

5 Dissassembly


Once your program is complete, remove the wrap from your injury site and disassemble the Pump and Therapy Units. Clean the device components and store away

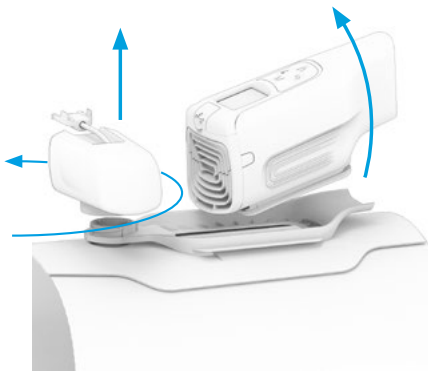

ProMOTION APP

1 Connecting


Turn the Therapy Unit on, and in the app, search for your device. Select your device from the list and wait for it to connect.

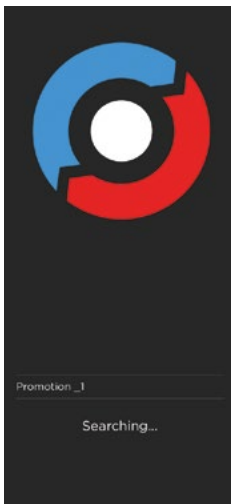

ProMOTION APP

2 Personal Information

Complete the Personal Information section as required.

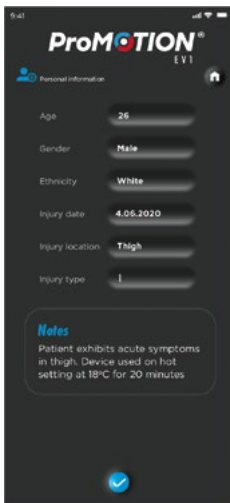

The screenshot shows the ProMOTION EV1 app interface. At the top, the status bar displays the time 9:41, signal strength, Wi-Fi, and battery icons. The app title "ProMOTION EV1" is centered at the top. Below the title, there is a "Personal information" section header with a person icon on the left and a home icon on the right. The form contains the following fields:


Age	26
Gender	Male
Ethnicity	White
Injury date	4.06.2020
Injury location	Thigh
Injury type	I

Below the form is a "Notes" section with a blue header. The text in the notes reads: "Patient exhibits acute symptoms in thigh. Device used on hot setting at 18°C for 20 minutes". At the bottom center of the screen is a blue circular button with a white checkmark, indicating that the form is complete.

ProMOTION APP

3 Program selection

Select 'Default Program' to choose one of our predefined programs or select 'New Program' to define your own.

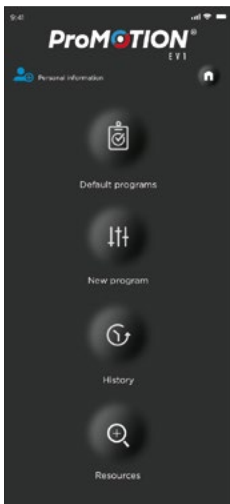

ProMOTION APP

4 Input

Enter the desired parameters for your program. Temperature, duration and compression can be adjusted to suit. Initiate the program once you are ready to go.



ProMOTION[®]
E V I

Customer Services

0845 901 0227

customerservices@swellaway.com

www.promotion.fitness